

Type 1036 In-Line Sanitary Pressure Gauge with Type 1037 Sanitary Instrument Fitting





FEATURES

Type 1036 Sanitary Gauge

- · All-welded stainless steel wetted surfaces
- Field liquid-fillable gauge case
- True Zero[™] pointer indication
- PowerFlex™ movement for extended life
- Easv Zero[™] external pointer adjustment standard
- Retrofits Anderson Instrument CPM design

Type 1037 Fitting

- Pipe O.D. sizes from 1/2" thru 2"
- 316L SS
- Electropolished 12-20RA internal surface finish
- Heat number of wetted materials stamped on each fitting


Designed specifically for the pharmaceutical, biotechnology and food processing markets, the Ashcroft® Type 1036 in-line sanitary pressure gauge and Type 1037 sanitary instrument fitting virtually eliminate process dead-leg. The design of the Type 1036 sanitary gauge and instrument fitting allows for the diaphragm of the gauge to be positioned at the gauge instrument fitting, eliminating the pocket or deadleg that may cause contamination.

The Type 1036 sanitary gauge and Type 1037 instrument fitting utilize a 1½" Tri-Clover-type mating connection. This feature offers flexibility to use the Ashcroft Type 1036 with the Type 1037 instrument fitting on sanitary instrument fittings for pipe sizes from 1/2" thru 2".

The Ashcroft Type 1036 sanitary pressure gauge can be auto-claved or sterilized(1), and cleaned (CIP) or steamed in place (SIP). The gauge comes standard with True Zero pointer indication. Since the Type 1036 has no dial stop pin, the gauge pointer always indicates a true pointer reading.

(1) Optional polysulfone window required.

HOW TO ORDER (Typical example)

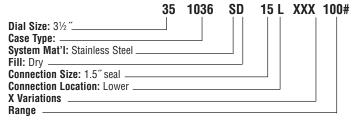

75-1037

HOW TO ORDER

© Ashcroft Inc. 2007

50 - 1/2" Tri-Clamp connection

75 – ¾" Tri-Clamp connection

10 – 1" Tri-Clamp connection

15 - 11/2"Tri-Clamp connection

20 – 2" Tri-Clamp connection





Type 1036 In-Line Sanitary Pressure Gauge with Type 1037 Sanitary Instrument Fitting

SPECIFICATIONS - Type 1036

Dial Size: 3½

Accuracy: 1.5% F.S. for pressure ranges 100 psi and above. 2.0% F.S. for vacuum, compound and ranges

below 100 psi. Add 0.5% if liquid-filled.

Case & Ring Material: 300 series stainless steel

Ring Type: Bayonet, removable

Bourdon Tube and Socket Material:
316 stainless steel

Diaphragm Material & Surface Finish:

316 stainless steel, 12 to 20RA

Diaphragm O-Ring: Buna-N(1) (standard)

Viton^{® (1)} (XVD), Teflon⁽¹⁾ (XTD)

Connection Location: Lower

System Fill: USP food grade glycerin Window Material: Polycarbonate

Pointer: Black-painted aluminum with external zero

adjust (Easy Zero™)

Dial: white with black markings **Movement:** 300 series stainless steel

Agency Approvals: 3A compliance to sanitary

standard 74-02

Ranges: 15# thru 600#, including compound

& vacuum

Clean or Steam in Place: (CIP or SIP): Temperature limit of 300°F (149°C)

Autoclave or Sterilize: Temperature limit of

300°F (149°C)

We recommend a polysulfone window for autoclave/sterilization. Specify the XPS variation.


SPECIFICATIONS - Type 1037	CODE
316L stainless steel construction	std.
Wetted parts electropolished to 12 to 20RA surface finish	std.
Heat number stamped on fitting	std.
Sizes:	
½" Tri-Clamp connection	50
¾" Tri-Clamp connection	75
1" Tri-Clamp connection	10
1½ Tri-Clamp connection	15

2" Tri-Clamp connection(1) To Ensure Cleanliness:


Prior to reinstalling the Type 1036 in to the Type 1037 instrument fitting, we recommend replacing the O-ring (P/N 185A106-75) Buna-N

20

185A106-75 Viton 185A106-77 Teflon

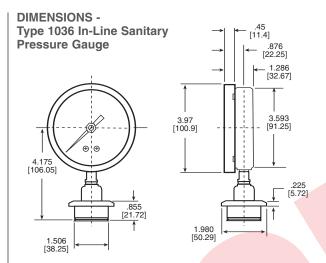

DIMENSIONS -Type 1037 Instrument Fitting

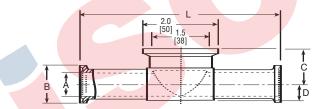

Inch/mm	Α	В	С	D	L
1/2"	.38	.98	.81	.31	5.25
50	(9.6)	(24.9)	(20.6)	(7.9)	(133.4)
3/4"	.63	.98	.93	.43	5.25
75	(15.9)	(24.9)	(23.6)	(10.9)	(133.4)
1″	.87	1.98	1.05	.55	4.75
10	(22.1)	(50.4)	(26.7)	(14.1)	(120.7)
11/2"	1.37	1.98	1.31	.81	5.50
15	(34.7)	(50.4)	(33.3)	(20.6)	(139.7)
2″	1.86	2.51	1.56	(1.06	(7.00
20	(47.2)	(62.8)	(39.6)	26.9)	177.8)